


Updating the EXA/J-Screen software


1. Unplug the USB flash drive from your EXA/J-Screen and plug into a computer connected to the web.
2. Create a backup copy of the content of your EXA/J-Screen USB flash drive on your computer.
3. Delete the content of your EXA/J-Screen USB flash drive.
4. Download the zip file from: <http://www.gemmoraman.com/Portals/0/Documents/EXA.zip>
5. Extract the content of the EXA.zip file in an empty folder on your PC. This file may show up automatically in your Downloads folder.
5. Copy the content of the folder to your EXA/J-Screen USB flash drive. If you previously saved your own spectra on the old USB flash drive, copy them into the "MySpectra" folder of the new USB flash drive. Create the MySpectra- folder manually to the root of USB flash drive in case it does not already exist.

Before starting the software updating on your EXA/J-Screen it's necessary to plug a mouse into the USB port of the unit to make the installation process easier.


6. Plug the USB flash drive into your EXA/J-Screen.
7. Turn on your EXA/J-Screen.

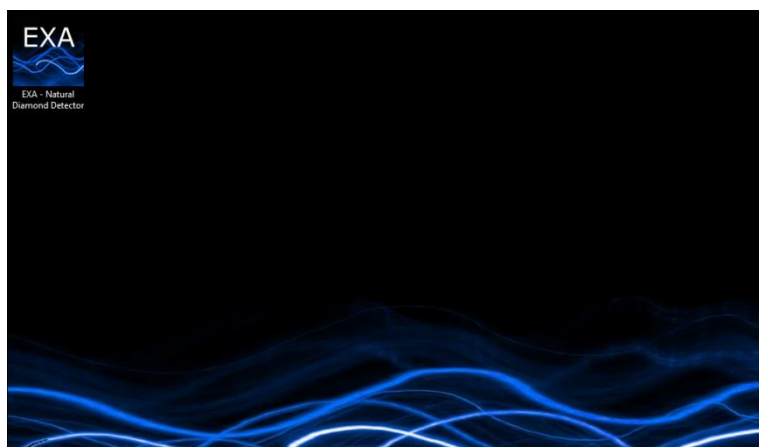


8. Tap the windows icon on the bottom-left corner of the screen

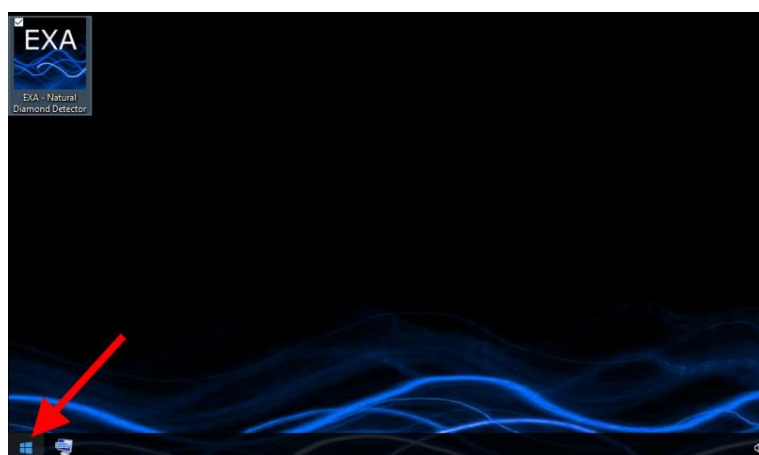


9. Open the Windows system submenu.
10. Launch the file explorer under Windows system.

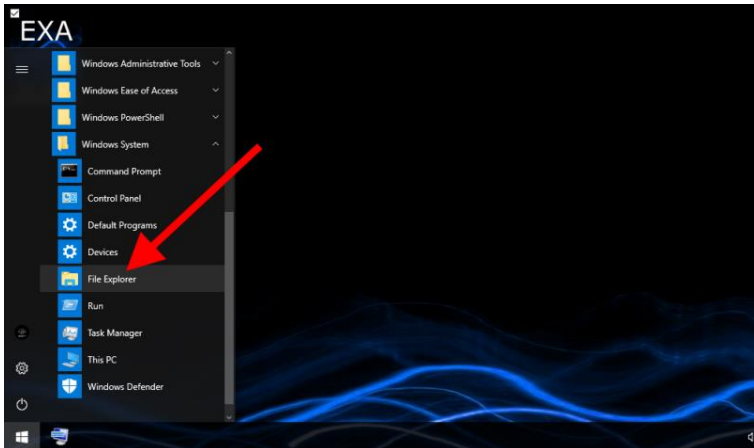


11. Open the content of the USB flash drive (D:)

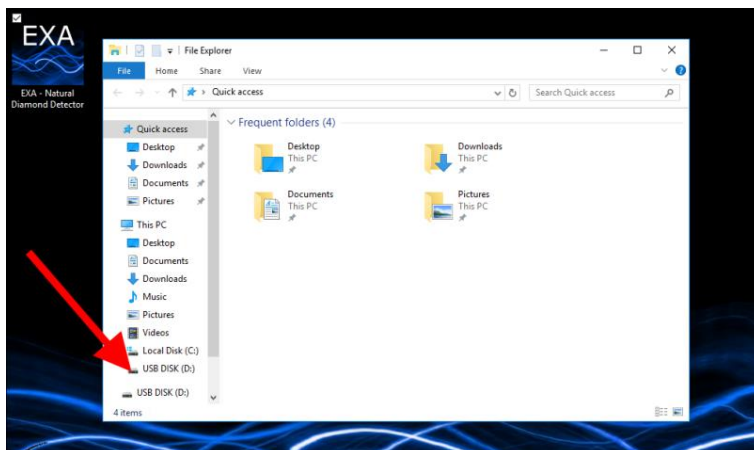

12. Launch the “EXA Upgrade.exe” file.

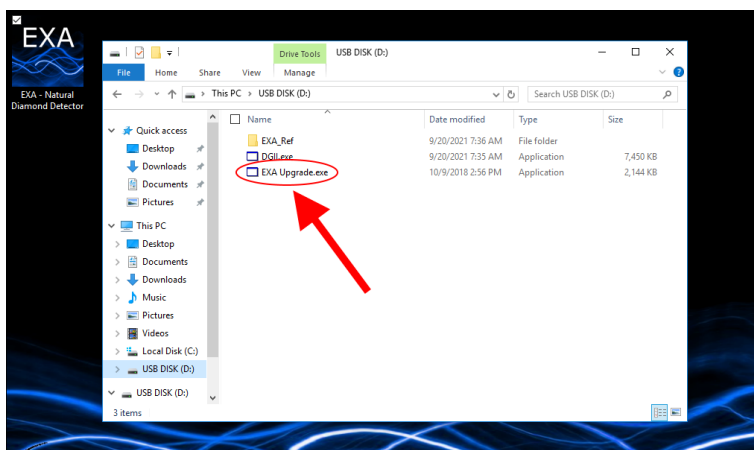

13. Click Start on the pop up menu.

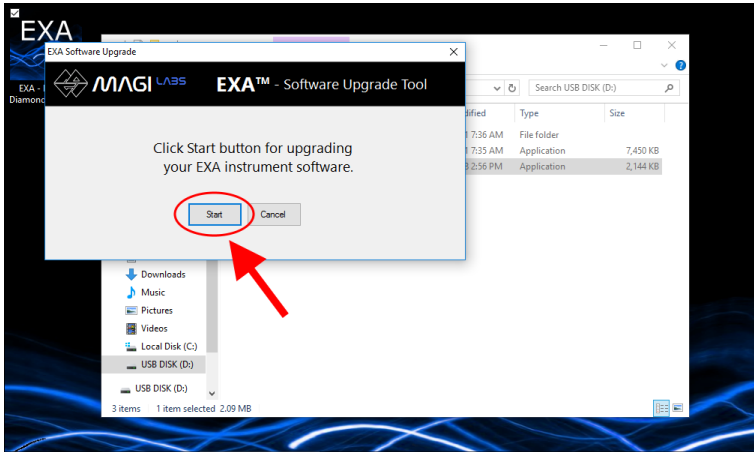

In Few seconds the software will be updated. Now you can start using the unit again by clicking on the EXA/J-Screen icon on the desktop.